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MOZOGRAU®	Technical File. Section 2: INFORMATION TO BE SUPPLIED BY THE MANUFACTURER	CE MARKING CE
	2.2.1 IFU (IMPLANT MOUNTS / DRIVERS)	Rev.: 01
		Date: 24/01/2025
		Page 4 of 10

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MOZOGRAU®	Technical File. Section 2: INFORMATION TO BE SUPPLIED BY THE MANUFACTURER	CE MARKING CE
	2.2.1 IFU (IMPLANT MOUNTS / DRIVERS)	Rev.: 01
		Date: 24/01/2025
		Page 5 of 10

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MOZOGRAU®	Technical File. Section 2: INFORMATION TO BE SUPPLIED BY THE MANUFACTURER	CE MARKING CE
	2.2.1 IFU (IMPLANT MOUNTS / DRIVERS)	Rev.: 01
		Date: 24/01/2025
		Page 6 of 10

ENGLISH

SCOPE

This instruction for use has its scope in implant mounts and drivers manufactured by Mozo Grau.

The implant mounts and drivers of TICARE implant systems are tools made of titanium or stainless steel intended to be used in the process of placing a Ticare Dental Implants for transfer the implant to the patient's mouth and its subsequent placement.


INDICATIONS FOR USE

TICARE Implant Mounts and drivers are intended to be used in the process of placing a TICARE Dental Implant to transfer the implant to the patient's mouth and its subsequent placement without compromising the sterility of the implant material

Function of use of implant mounts and drivers: **Ensure correct transport of the implant within the vial and subsequent placement in the patient's mouth.**

Characteristics:

- They are made of stainless steel and TiAl6V4 ELI (Grade 5 Ti) Ergitan 3.7165 MG
- Material testing in accordance with ASTM A582, ASTM F136 and ASTM F899.
- Color coding, indicating the diameter of the implant being transported:

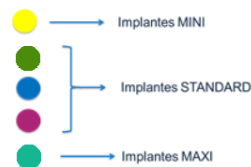USER

Implant mounts and drivers are designed to be used by qualified dental surgeons or dental health professionals who have received additional training in dental implant procedures (dental implantology).

INTENDED TARGET GROUPS

Implant mounts and drivers are intended to be used in edentulous and partially dentate patients, with consolidated bone and completed growth period, and who are subject to dental implant treatments or aesthetic oral rehabilitation.

CONTRAINDICATIONS

Dental implant procedures are contraindicated for patients who are not medically fit to undergo an oral surgery procedure.

Contraindications include:

1. Patients allergic or hypersensitive to the raw materials used.
2. Patients with mental problems that prevent patient cooperation.
3. Patients who abuse drugs and alcohol.
4. Patients with uncontrolled diabetes or blood disorders
5. Dental implant procedures are not recommended in pediatric patients until the growth of the maxillary bone is complete.

MOZOGRAU®	Technical File. Section 2: INFORMATION TO BE SUPPLIED BY THE MANUFACTURER	CE MARKING CE
	2.2.1 IFU (IMPLANT MOUNTS / DRIVERS)	Rev.: 01
		Date: 24/01/2025
		Page 7 of 10

Allergen information:

Materials used for the manufacture of implant mounts and drivers:

- (TiAl6V4 ELI (Grade 5 Ti) Ergitan 3.7165 MG)
- AISI 303 (Stainless steel)
- AISI 440B (Martensitic stainless steel)

WARNINGS

Do not use if individual packaging is open or damaged.

Products must be protected against aspiration when handled intraorally. Aspiration of products may cause unintended physical injury or infection. Always inspect instruments before use.

Instruments must be used only by trained and specialized dental personnel. For safe and effective use of the system, it is strongly suggested that specific training be undertaken, as well as the study of appropriate books. THESE INSTRUCTIONS ARE NOT INTENDED TO BE A SUBSTITUTE FOR PROPER TRAINING.

Instruments must be handled carefully, avoiding the use of damaged products as this may compromise their effectiveness.

After each use, clean and disinfect the product. For proper maintenance of the instruments, the cleaning and sterilization instructions must be followed. Do not store wet instruments after the cleaning and/or sterilization process.

TICARE implant mounts do not require special storage conditions, as the raw material is stable under atmospheric pressure and temperature conditions. Transport and storage temperature <50°C and humidity <50%. TICARE implant mounts and drivers are stable under atmospheric transport and storage conditions.

For proper fit and positioning of TICARE implant mounts, use only original TICARE implants.

PRECAUTIONS

The product is supplied NON-STERILE and must be used under aseptic conditions.

The products included in this instruction can be reused. Reuse of these products is possible, provided that strict control is exercised over their cleaning and sterilization.

Under normal conditions, TICARE implant mount and drivers have an established useful life of 10 years.

Special care is recommended when treating patients who have local or systemic risk factors that may interfere with the healing process of bone or soft tissue. These risk factors include, among others, smoking, poor oral hygiene, uncontrolled diabetes, pregnancy, oral radiotherapy, infections in the neighboring bone and certain therapies, including steroids and chemotherapy.

Advise patients who smoke to stop smoking before undergoing dental implant procedures.

Dental implant procedures are associated with ancillary medical procedures (e.g. x-rays, administration of anesthetics, antibiotic treatments, etc.) which may be contraindicated in some medical conditions such as pregnancy or breastfeeding, among others. It is the final responsibility of the dentist or qualified professional to judge whether patients are medically fit to undergo an oral surgery procedure. Regular dental check-ups and maintaining good oral hygiene are crucial to minimize adverse effects.

Take special care when transporting the implant, as there is a risk of it becoming dislodged and falling out. If the O-ring is in poor condition, there is an increased risk of the implant falling out during transport. If the implant falls out, do not use it and replace it with a new one. Due to the small size of the devices, care must be taken to ensure that the patient does not swallow or aspirate them. It is appropriate to use specific support tools to prevent aspiration of loose parts (e.g. gauze, dental dam or throat guard).

MOZOGRAU®	Technical File. Section 2: INFORMATION TO BE SUPPLIED BY THE MANUFACTURER	CE MARKING CE
	2.2.1 IFU (IMPLANT MOUNTS / DRIVERS)	Rev.: 01
		Date: 24/01/2025
		Page 8 of 10

CLEANING AND STERILIZATION

A) Treatment instructions/ warnings

To avoid damage, do not use metal brushes or cleaning sponges.

Only use cleaning and disinfecting solutions with a pH value between 4.5 and 10. These instruments are supplied non-sterile and must be cleaned and sterilized before use.

B) Post processing restriction

The end of the product's useful life is normally marked by wear and tear resulting from incorrect use and handling.

C) Cleaning, maintenance and sterilization of instruments

Surgical instruments can deteriorate if proper maintenance is not carried out. At Mozo Grau we recommend our customers carry out appropriate cleaning and maintenance procedures before and after use to ensure a good useful life of the instrument.

C.1) During the use of surgical material: dirty products must be separated so as not to contaminate the instrument tray that is being used. Dirty products should be cleaned and replaced in the instrument tray before proceeding to sterilize it with the necessary components.

C.2) Treatment after surgical use:

Major dirt should be removed from the products immediately after use. To do this, use pure water or a solution with disinfectant, which should not contain aldehyde (otherwise the blood would settle into the product).

1 Selection of the product for cleaning/disinfection. It would be recommended to use a neutral/enzymatic cleaner (pH allowed between 6 and 10): It should not contain bleach, organic solvents (Alcohols, Ethers, etc.), oxidants (e.g., hydrogen peroxide), or halogens (chlorine, bromine, iodine). Be careful with tap water that is too chlorinated. Avoid solution with aldehyde, due to its ability to bind blood to the surgical material. No specific cleaning or disinfecting agent is recommended, but it is necessary to observe the instructions on concentrations and immersion times of the devices in the cleaning or disinfecting solutions indicated by the manufacturer.

2 Selection of material to carry out cleaning. For manual dirt removal, use only nylon brushes; Never use metal brushes or steel wool.

3 How to carry out cleaning/disinfection. Carefully remove all postoperative residues (Blood, bone...). Scrub the devices with a soft-bristled nylon brush in tap water until all dirt has been removed. Immerse the devices in an ultrasonic bath with a cleaning solution at the concentration and temperature specified in the detergent manufacturer's instructions.

4 Rinse for one minute, preferably with purified or distilled water until the cleaning solution has been removed. If there are still traces of dirt after cleaning, repeat the same procedure.

Then, immerse in a bath with disinfectant solution at the concentration and temperatures indicated by the manufacturer. Make sure all ducts are filled with disinfectant solution and all surfaces are covered with the solution.

5 Rinse with purified or distilled water until all the disinfectant solution has been removed.

6 Dry the devices immediately after the cleaning and disinfection process with single-use wipes. Never store instruments while they are still wet.

MOZOGRAU®	Technical File. Section 2: INFORMATION TO BE SUPPLIED BY THE MANUFACTURER	CE MARKING CE
	2.2.1 IFU (IMPLANT MOUNTS / DRIVERS)	Rev.: 01
		Date: 24/01/2025
		Page 9 of 10

D) How to perform sterilization

Moist heat sterilization should be performed in an autoclave at 134° (273°F) for 5 minutes. Make sure that the elements inside the autoclave are not rusty. Surgical material must be sterilized inside sterilization bags. Do not sterilize materials of different types in contact. The exposure temperature of the products and sterilization trays will not exceed 137°C.

Only for the United States: The validated parameters are T^a 270°F / 132°C for 4 minutes and a drying time of 20 minutes.

PACKAGING AND STORAGE

The product is delivered clean, but not sterile. All information about the product is on the label. Product packaging is done individually. No special storage conditions are needed since the raw materials used are stable under normal pressure and temperature conditions. It is recommended to store between temperatures between 18°C and 50°C and a humidity limit below 50%. Store the set at a moderate temperature in a dry, dust-free and ventilated place where corrosive vapors are not produced.

DISPOSAL

The product should be disposed of in an environmentally friendly manner in accordance with local regulations. Hazardous waste of contaminated devices or sharp objects should be disposed of in suitable containers which meet the specific technical requirements.

DETECTION AND NOTIFICATION OF INCIDENTS WITH HEALTH PRODUCTS

Important notice: Any patient or healthcare personnel who detects any serious incident related to the Ticare dental implant must notify the manufacturer, Mozo Grau S.A. and the competent authority of the Member State in which the user and/or patient are established.

MOZO GRAU®	Technical File. Section 2: INFORMATION TO BE SUPPLIED BY THE MANUFACTURER	CE MARKING CE
	2.2.1 IFU (IMPLANT MOUNTS / DRIVERS)	Rev.: 01
		Date: 24/01/2025
		Page 10 of 10

LABELLING SYMBOLS

SYMBOLS	DESCRIPTION
	Reference of the producto
	Batch code
	Manufacture date
	Manufacturer
	Medical device
	Non sterile
	CE Mark
	Consult electronic instructions for use
	Unique device identifier

MANUFACTURE

Mozo Grau S.A.
 C/ Santiago López González, N°7
 47197 Valladolid (Spain)
 Tel: +34 983 309 602
 e-mail: info@ticareimplants.com
www.ticareimplants.com