MOZOGRAU®	Technical File. Section 2: INFORMATION TO BE SUPPLIED BY THE MANUFACTURER	CE MARKING CE
	2.2.4 IFU (PROFILE DRILLS/ FRESAS PERFILADORAS)	Rev.: 03
		Date: 18/12/2024
		Página 1 de 17

ESPAÑOL

ALCANCE

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ENGLISH

SCOPE

This instruction for use has its scope in Dental profile drills and countersink manufactured by Mozo Grau.

The Dental profile Drills and Countersink of TICARE implant systems are mechanical cutting tools, made of stainless steel. They are helical in shape and are used in conjunction with a handpiece connected to a motor.

The part connecting the contra-angle is manufactured in accordance with the current ISO 1797.

All rotatory and oscillating instruments manufactured by Mozo Grau are according to the current ISO 3823 standard


INDICATIONS FOR USE

These dental instruments are indicated for the drilling and / or preparation of the bone in oral implantology processes for the placement of a dental implant.

Function of use of the Dental Drills: Preparation of the socket to place the implant. To clean the bone of the crestal area of the implant and facilitate the insertion of the prosthesis.

Characteristics:

- They are stainless steel
- Laser marking indicates the insertion limit of the drill bit.
- **Profile drills:** To widen the implant bed in the cortical area to accommodate the micro-thread area without overcompressing the bone, in Inhex implants (Internal Connection).
 - Color coding, which indicates the diameter of the implant to be placed and laser marking on the stem to identify the diameter:

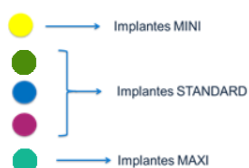- **Countersink:** Only for Ossoeous implants (External Connection)
 - Only for 3,4 y 3,75 (Diameter of 4.1 mm)

USER

Profile drills and Countersink are designed to be used by qualified dental surgeons or dental health professionals who have received additional training in dental implant procedures (dental implantology).

INTENDED TARGET GROUPS

Profile drills and Countersink are intended to be used in edentulous and partially dentate patients, with consolidated bone and completed growth period, and who are subject to dental implant treatments or aesthetic oral rehabilitation.

CONTRAINDICATIONS

Lack of quality and quantity of remaining bone, infection, and local or systemic factors that may interfere with the healing process of bone or soft tissue (poor oral hygiene, uncontrolled diabetes, pregnancy, oral radiotherapy, i.a.) are some potential causes of stability loss.

Dental implant procedures are contraindicated for patients who are medically unfit to undergo an oral surgical procedure

MOZOGRAU®	Technical File. Section 2: INFORMATION TO BE SUPPLIED BY THE MANUFACTURER	CE MARKING CE
	2.2.4 IFU (PROFILE DRILLS/ FRESAS PERFILADORAS)	Rev.: 03
		Date: 18/12/2024
		Página 11 de 17

Contraindications include:

1. Patients allergic or hypersensitive to the raw materials used. Stainless steel
2. Patients with mental problems that prevent the patient's cooperation.
3. Patients who abuse drugs and alcohol.
4. Patients with uncontrolled diabetes or blood disorders
5. Implant procedures with dental drills in pediatric patients is not recommended until jaw bone growth is completed.

Allergen information:

The materials used for the manufacture of dental surgical burs:

- 465 KL (Stainless Steel) (C, Si, Mn, P, S, Cr, Mo, Ni, Ti, Fe)

Products that have changed their properties with use, such that they are not safe for the functions indicated, will not be used.

WARNINGS

Do not use if the individual package is open or damaged.

Products must be secured against aspiration when handled intraorally. Aspiration of products may lead to infection or unplanned physical injury. Always inspect instruments before use. Avoid the mandibular nerve canal during implant bed preparation and implant insertion. Nerve damage may result in anesthesia, paresthesia, or dysesthesia. Do not exceed the recommended insertion torques as this might cause bone necrosis.

The instruments used exclusively by personnel trained and specialized in dentistry. For the safety and effectiveness of the use of the system, it is strongly suggested that specific training be carried out, as well as the study of appropriate books. THESE INSTRUCTIONS ARE NOT INTENDED TO BE A SUBSTITUTE FOR ADEQUATE TRAINING.

Instrument handling must be carried out carefully, avoiding the use of damaged products as this may compromise their effectiveness. Check the condition of the cutting edges. If there is any sign of wear, discard the surgical drills.

After each use, clean and disinfect the product. For good maintenance of the instruments, the cleaning and sterilization instructions must be followed. Do not store wet instruments after the cleaning and/or sterilization process.

Failure to use radiographic measurements when planning surgery and the drilling sequence can result in permanent damage to vital oral structures

Drilling beyond the depth intended for lower jaw surgery may result in numbness to the lower lip or chin or lead to hemorrhage in the floor of the mouth. There should be at least 2mm of bone under the actual implant. For example, to place an 8mm-length implant, 10 mm of bone height is required. The crest profile should be carefully examined to calculate the insertion angle to achieve parallelism with additional implants and natural dentition.

Failure to refer to anatomic knowledge and specific pre-operative radiographs may result in damage to vital oral structures when drilling in the jaw bone.

PRECAUTIONS

All product is supplied NON-STERILE, and must be used under aseptic conditions.

The products included in this instruction can be reused. The reuse of these products is possible, as long as strict control is maintained over their cleaning and sterilization.

It is recommended that after 20 uses and their corresponding sterilizations after each surgery, the uncoated stainless steel instruments be replaced to guarantee the cutting edge.

MOZOGRAU®	Technical File. Section 2: INFORMATION TO BE SUPPLIED BY THE MANUFACTURER	CE MARKING CE
	2.2.4 IFU (PROFILE DRILLS/ FRESAS PERFILADORAS)	Rev.: 03
		Date: 18/12/2024
		Página 12 de 17

Stainless steel is a resistant material, although it is not rust-free if it is not properly maintained.

Careful clinical and radiological examination of the patient is essential to treatment planning prior to surgery.

Special care is advised when treating patients who have local or systemic risk factors that may interfere with the healing process of bone or soft tissue. These risk factors include, i.a., smoking, poor oral hygiene, uncontrolled diabetes, pregnancy, oral radiotherapy, infections in neighbouring bone and certain therapies including steroids and chemotherapy. Bruxism or other jaw disorder may impede the success of dental implant procedures and may require the consideration of alternative treatment options.

Advise patients who smoke to stop smoking before undergoing dental implant procedures.

Dental implant procedures are associated with auxiliar medical acts (ex: x-rays, administration of anaesthetics, antibiotic treatments, etc.) that may be contra-indicated in some medical conditions such as pregnancy or breastfeeding amongst others. It is the final responsibility of the dental practitioner or qualified professional to judge whether the patients are medically fit or unfit to undergo an oral surgical procedure.

Regular dental check-ups and maintaining good oral hygiene are crucial to minimize any adverse effects..

CLEANING AND STERILIZATION

A) Treatment instructions/warnings

To avoid damage, do not use metal brushes or cleaning sponges.

Only use cleaning and disinfecting solutions with a pH value between 4.5 and 10. These instruments are supplied non-sterile and must be cleaned and sterilized before use.

B) Post processing restriction

The end of the product's useful life is normally marked by wear and tear resulting from incorrect use and handling.

C) Cleaning, maintenance and sterilization of instruments

Surgical instruments can deteriorate if proper maintenance is not carried out. At Mozo Grau we recommend our customers carry out appropriate cleaning and maintenance procedures before and after use to ensure a good useful life of the instrument.

C.1) During the use of surgical material: dirty products must be separated so as not to contaminate the instrument tray that is being used. Dirty products should be cleaned and replaced in the instrument tray before proceeding to sterilize it with the necessary components.

C.2) Treatment after surgical use:

Major dirt should be removed from the products immediately after use. To do this, use pure water or a solution with disinfectant, which should not contain aldehyde (otherwise the blood would settle into the product).

1 Selection of the product for cleaning/disinfection. It would be recommended to use a neutral/enzymatic cleaner (pH allowed between 6 and 10): It should not contain bleach, organic solvents (Alcohols, Ethers, etc.), oxidants (e.g., hydrogen peroxide), or halogens (chlorine, bromine, iodine). Be careful with tap water that is too chlorinated. Avoid solution with aldehyde, due to its ability to bind blood to the surgical material. No specific cleaning or disinfecting agent is recommended, but it is necessary to observe the instructions on concentrations and immersion times of the devices in the cleaning or disinfecting solutions indicated by the manufacturer.

MOZOGRAU®	Technical File. Section 2: INFORMATION TO BE SUPPLIED BY THE MANUFACTURER	CE MARKING CE
	2.2.4 IFU (PROFILE DRILLS/ FRESAS PERFILADORAS)	Rev.: 03
		Date: 18/12/2024
		Página 13 de 17

2 Selection of material to carry out cleaning. For manual dirt removal, use only nylon brushes; Never use metal brushes or steel wool.

3 How to carry out cleaning/disinfection. Carefully remove all postoperative residues (Blood, bone...), in the case of products with cavities, such as drills with irrigation, rinse all cavities of the products using a disposable syringe. Scrub the devices with a soft-bristled nylon brush in tap water until all dirt has been removed. Immerse the devices in an ultrasonic bath with a cleaning solution at the concentration and temperature specified in the detergent manufacturer's instructions.

4 Rinse for one minute, preferably with purified or distilled water until the cleaning solution has been removed. If there are still traces of dirt after cleaning, repeat the same procedure.

Then, immerse in a bath with disinfectant solution at the concentration and temperatures indicated by the manufacturer. Make sure all ducts are filled with disinfectant solution and all surfaces are covered with the solution.

5 Rinse with purified or distilled water until all the disinfectant solution has been removed.

6 Dry the devices immediately after the cleaning and disinfection process with single-use wipes. Never store instruments while they are still wet.

D) How to perform sterilization

Moist heat sterilization should be performed in an autoclave at 134° (273°F) for 5 minutes. Make sure that the elements inside the autoclave are not rusty. Surgical material must be sterilized inside sterilization bags. Do not sterilize materials of different types in contact. The exposure temperature of the products and sterilization trays will not exceed 137°C.


Only for the United States: The validated parameters are T^a 270°F / 132°C for 4 minutes and a drying time of 20 minutes.

SURGICAL PROCEDURES


1) Inhex and Inhex hybrid milling protocol

MOZOGRAU®	Technical File. Section 2: INFORMATION TO BE SUPPLIED BY THE MANUFACTURER	CE MARKING CE
	2.2.4 IFU (PROFILE DRILLS/ FRESAS PERFILADORAS)	Rev.: 03
		Date: 18/12/2024
		Página 14 de 17

STANDARD DRILLING SEQUENCE

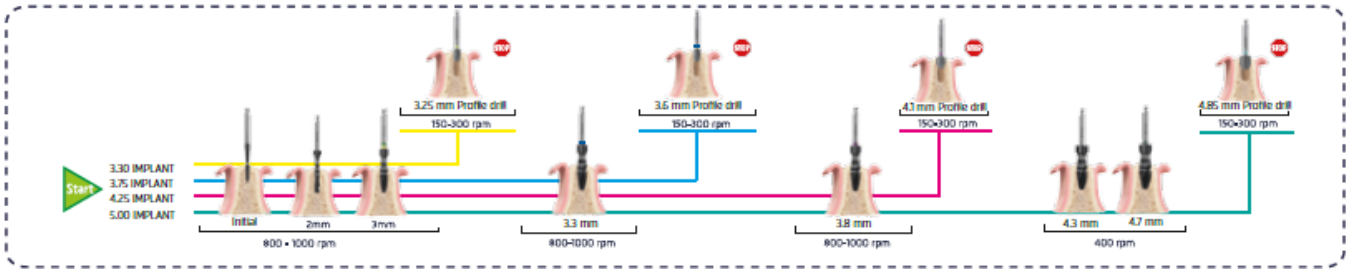
HARD BONE DRILLING SEQUENCE

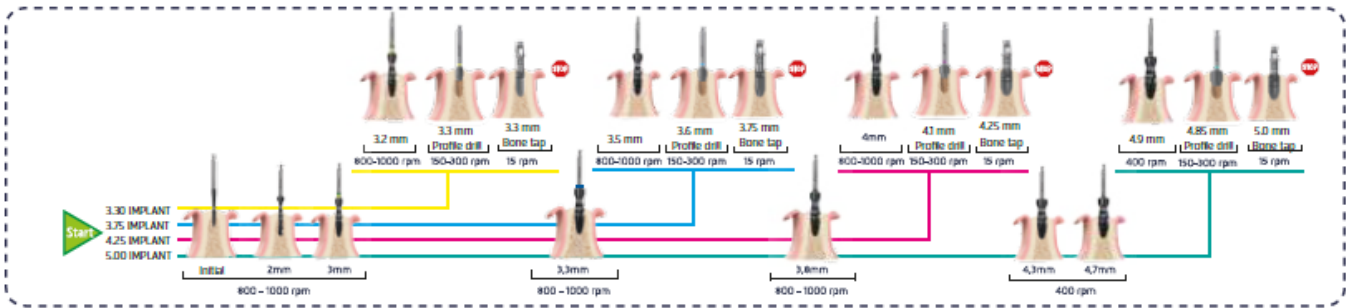
Profile drills should not be used for 6 mm long implants

2) Inhex Quattro and Inhex Quattro hybrid milling protocol

MOZOGRAU®	Technical File. Section 2: INFORMATION TO BE SUPPLIED BY THE MANUFACTURER	CE MARKING CE
	2.2.4 IFU (PROFILE DRILLS/ FRESAS PERFILADORAS)	Rev.: 03
		Date: 18/12/2024
		Página 15 de 17

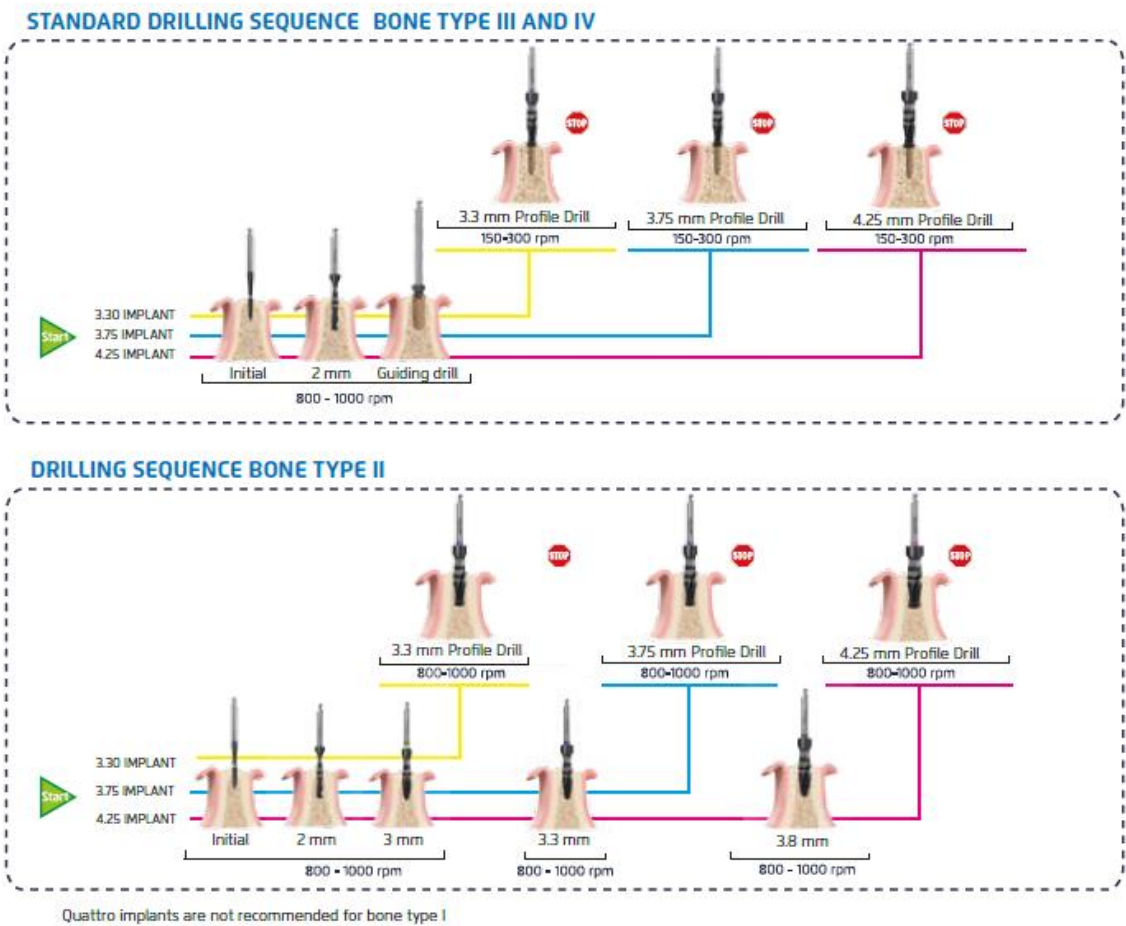
3) Osseous milling protocol

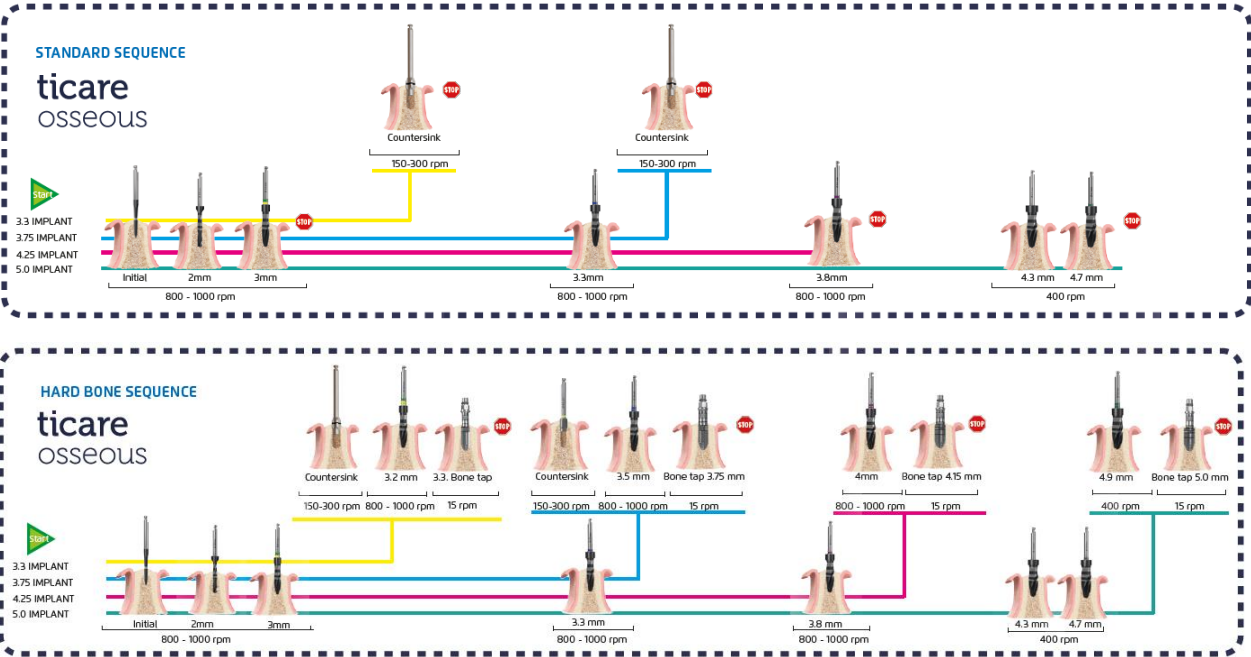MOZOGRAU®	Technical File. Section 2: INFORMATION TO BE SUPPLIED BY THE MANUFACTURER	CE MARKING CE
	2.2.4 IFU (PROFILE DRILLS/ FRESAS PERFILADORAS)	Rev.: 03
		Date: 18/12/2024
		Página 16 de 17

4) Osseous Quattro milling protocol

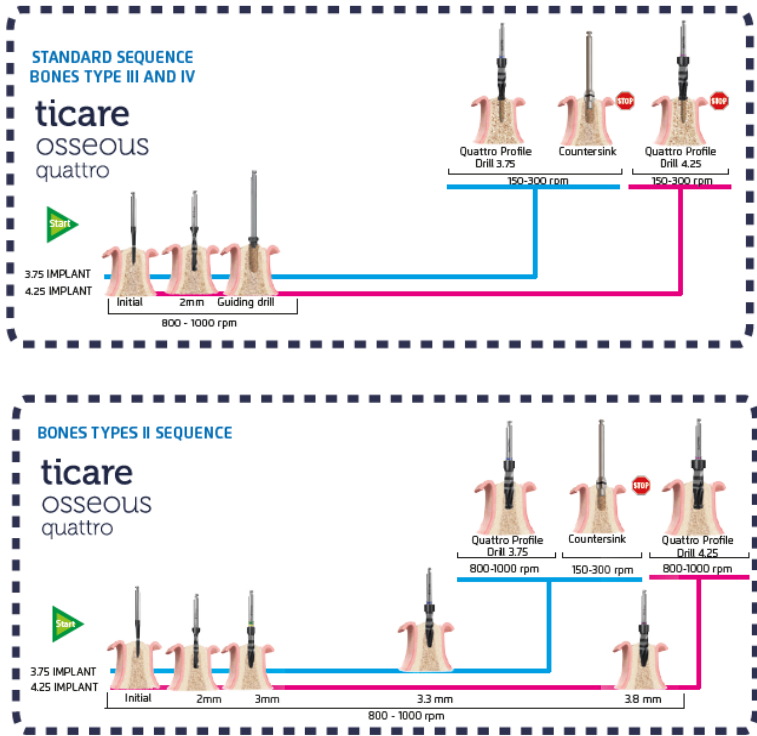PACKAGING AND STORAGE

The product is delivered clean, but not sterile. All information about the product is on the label. Product packaging is done individually. No special storage conditions are needed since the raw materials used are stable under normal pressure and temperature conditions. It is recommended to store between temperatures between 18°C and 50°C and a humidity limit below 50%. Store the set at a moderate temperature in a dry, dust-free and ventilated place where corrosive vapors are not produced.

DISPOSAL

The product should be disposed of in an environmentally friendly manner in accordance with local regulations. Hazardous waste of contaminated devices or sharp objects should be disposed of in suitable containers which meet the specific technical requirements.

DETECTION AND NOTIFICATION OF INCIDENTS WITH HEALTH PRODUCTS

Important notice: Any patient or healthcare personnel who detects any serious incident related to the Ticare dental implant must notify the manufacturer, Mozo Grau S.A. and the competent authority of the Member State in which the user and/or patient are established.

LABELLING SYMBOLS

SYMBOLS	DESCRIPTION
REF	Reference of the producto
LOT	Batch code
	Manufacture date

MOZOGRAU®	Technical File. Section 2: INFORMATION TO BE SUPPLIED BY THE MANUFACTURER	CE MARKING CE
	2.2.4 IFU (PROFILE DRILLS/ FRESAS PERFILADORAS)	Rev.: 03
		Date: 18/12/2024
		Página 17 de 17

	Manufacturer
	Medical device
	Non sterile
	CE Mark
	Consult electronic instructions for use
	Unique device identifier

MANUFACTURE
Mozo Grau S.A.
C/ Santiago López González, N°7
47197 Valladolid (Spain)
Tel: +34 983 309 602
e-mail: info@ticareimplants.com
www.ticareimplants.com